




Sturdy and mess free bowl!



1



2

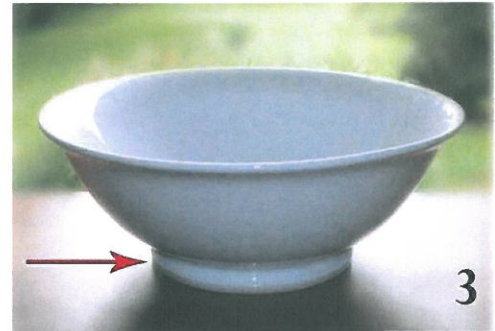

3

Materials:

Bowl with raised bottom to hold pennies (see picture 3)

Pennies

Plastic container top

Glue or tape

Directions:

1. Turn bowl upside down and place pennies into raised bottom of bowl (however many will fit is fine).
2. Apply glue or tape to raised edges of bottom of bowl. Also apply glue or tape to center of plastic container top (make sure edges of top are facing up). Place glued section of container top onto bottom of bowl so they are firmly connected. Apply more glue or tape as needed (see picture 2 for guidance).
3. Turn bowl over so that plastic container is on bottom and your device is ready to use (see picture 1)!

Idea submitted by: Pennsylvania's Initiative on Assistive Technology (PIAT)

For more information: ATinfo@temple.edu

<http://disabilities.temple.edu/programs/assistive/piat/>